



 PRODUCT CATALOGUE

2023

→ INDEX

6-7	PIEZOSURGERY® <i>touch</i>
8-9	PIEZOSURGERY® <i>white</i>
10	enzymec
11	universal cart
12-13	PIEZOSURGERY® spare parts
14	insert sets
15	basic set/osteotomy set
16	sinus lift lateral set/piezo lift set
17	sinus lift set/ implant prep set pro
18	extraction set/ explantation set
19	resective perio set/retro surgical set
20-21	indications
22	trays
23-25	implant prep inserts
26-29	osteotomy inserts
30-31	osteoplasty inserts
32	elevation inserts/periosteum inserts
33	endo retro inserts
34	extraction inserts/explantation inserts
35	periodontal inserts
36	implant cleaning inserts
37	crown prep tips

38-39 MULTIPIEZO

40-51 MICROPIEZO

42 COMPACT PIEZO LED/COMPACT PIEZO P2K

43-45 ultrasound accessories

46 insert sets

47 set maintenance/set scaling

48 set perio maintenance/set perio treatment

49 set implant cleaning/set endo revision

50 scaling inserts

51 perio universal inserts

52 perio anatomic inserts

53 implant cleaning inserts

54 endo inserts

55 endo revision inserts

56 endo retro inserts/ extraction inserts

57 crown prep inserts

→ **ULTRASOUND**

→ AIR-POLISHING

58-59 COMBI *touch*

60-61 STARJET

62-63 prophylaxis powders

64 universal cart

66-69 spare parts

→ LED-CURING

70-71 starlight uno

72-73 starlight pro

74-75 starlight s⁺/starlight sx5

76 spare parts

78-79 appropriate settings for PIEZOSURGERY® inserts

80-81 appropriate settings for ultrasound inserts

82-83 mectron online education

84-87 distributors

We reserve the right to make technical changes and changes in the scope of the supply. No liability can be accepted for any errors. The reference numbers of the units are meant for standard packagings; they may differ in some countries and according to different voltages.

PIEZOSURGERY® *touch*

The future of bone surgery? Get in touch.



PIEZOSURGERY® *touch*

Equipment set:

- 1 PIEZOSURGERY® *touch* unit
- 1 LED-handpiece, complete with cord
- 1 kit of 8 tubes for pump
- 1 torque wrench K8
- 1 foot pedal
- power cord and suitcase

Ref. no.

05120061



→ TECHNICAL DATA

sizes (l x w x h in mm):

300 x 250 x 95 mm

weight:

3.2 kg

operating frequency:

24.000 – 36.000 Hz

power types:

6 functions

IMPLANT

CORTICAL

CANCELLOUS

SPECIAL

PERIO

ENDO

peristaltic pump capacity:

adjustable on the touch screen:

ENDO / PERIO - 7 levels: from 0 to 6

(from 0 to approximately 75 ml/min)

SPECIAL / CANCELLOUS / CORTICAL /

IMPLANT - 6 capacity levels: from 1 to 6

(from 8 to approximately 75 ml/min)

power supply voltage:

digital 100-240 VAC, 50/60 Hz

max.power absorbed: 120 VA

warranty:

2 years, handpiece 1 year

PIEZOSURGERY® *white*

Equipment set:

1 PIEZOSURGERY® *white* unit

1 kit of 8 tubes for pump

1 torque wrench K8

1 foot pedal

power cord and suitcase
(handpiece and inserts not included)

Ref. no. 05120063

LED-HANDPIECE COMPLETE WITH CORD AND TRAY

- swivel-type LED light can be directed at the insert edge
- handpiece and handpiece cord (including the irrigation line) are fully sterilizable all at once

Ref. no. 03120134



HANDPIECE COMPLETE WITH CORD AND TRAY

handpiece and handpiece cord (including the irrigation line) are fully sterilizable all at once

Ref. no. 03120217



PIEZOSURGERY® *white*

Get started in bone surgery!





→ TECHNICAL DATA

sizes (l x w x h in mm):

300 x 250 x 95 mm

weight:

3.2 kg

operating frequency:

24.000 – 36.000 Hz

power types:

6 functions

IMPLANT

CORTICAL

CANCELLOUS

SPECIAL

PERIO

ENDO

peristaltic pump capacity:

adjustable on the touch screen:

ENDO / PERIO - 7 levels: from 0 to 6

(from 0 to approximately 75 ml/min)

SPECIAL / CANCELLOUS / CORTICAL /

IMPLANT - 6 capacity levels: from 1 to 6

(from 8 to approximately 75 ml/min)

power supply voltage:

digital 100-240 VAC, 50/60 Hz

max. power absorbed: 120 VA

warranty:

2 years, handpiece 1 year

→ ENZYMEC

Enzymatic detergent

ENZYMEC is intended for use in cleaning of dental and general surgical instruments, and can be used during pre-cleaning in automatic equipments, ultrasonic cleaners and washer disinfectors according to the instructions and recommendations of the equipment manufacturer.

→ ENZYMEC

Product	Ref. no.
Enzymec 4-bottle pack (4 liters)	03900002

- enzymatic solution for efficient removal of organic residues
- specifically dedicated to the “clean” function of all PIEZOSURGERY® devices
- with practical dosage cap



UNIVERSAL CART

For daily use

universal cart – maximum mobility for all your table-top devices. All supply lines are hidden in the easily accessible cable room. The power supply is brought at the tray's height, for an easy connection to your devices.



→ SIZES:

total (l x w x h in mm):

480 x 460 x 850 mm

upper tray:

480 x 370 in 770 mm height

middle tray:

420 x 310 in 440 mm height

lower tray:

480 x 370 in 65 mm height

weight:

29 kg

→ UNIVERSAL CART

Product

universal cart

Ref. no.

03540011



surgical tray



torque wrench K8



irrigation set



PIEZOSURGERY® insert box

SPARE PARTS

High quality spare parts for professionals



DPS adaptor



PIEZOSURGERY®
insert adaptor



ADAPTORS FOR THERMODISINFECTION

Product

Ref. no.

DPS Adaptor (for PIEZOSURGERY® dental irrigation line) *

04610006

spare filter for DPS Adaptor *

04590006

PIEZOSURGERY® insert Adaptor *

04610010

PIEZOSURGERY® COMMON SPARE PARTS

surgical tray

04610004

dynamometric torque wrench K8

02900109

dynamometric torque wrench K7 (for tipholder DB2)

02900081

irrigation set (50 pieces)

03230006

PIEZOSURGERY® insert box

02900111

insert tray

02900037

insert box with bone expander tray

02900107

insert box with implant insert tray

02900108

bag support rod

01380002

PIEZOSURGERY® *touch* AND PIEZOSURGERY® *white*

Product

LED-handpiece complete with cord (PS *touch* and PS *white*)

PIEZOSURGERY® *white* handpiece complete with cord

handpiece LED-cone

handpiece metal cone

handpiece holder

O-ring set for handpiece holder

mobile handpiece holder

kit of 8 tubes for pump

bone grafting kit (liquid flow rate reducer, equipped with 8 silicone tubes and 1 tube-tube connector)

protection foils

(pack of 10 pcs. - each in a single sterile packaging)

footswitch (PS *touch*)

footswitch (PS *white*)

PIEZOSURGERY® *touch* and *white* suitcase

PIEZOSURGERY® I + II + 3

handpiece complete with cord

handpiece holder (PSII + 3)

kit of 8 tubes for pump (PSI)

kit of 8 tubes for pump (PSII + 3,
for handpiece cords with metal connector)

kit of 8 tubes for pump (PSII + 3,
for handpiece cords with plastic connector)

peristaltic pump rolling (PSII + 3)

tube-tube connector

footswitch (PSI + II)

footswitch (PS3)

PIEZOSURGERY® 3 suitcase

Ref. no.

03120134

03120217

03020158

03020171

02900105

03020160

03110015

02900104

02900106

03230012

02900102

02900114

04440006

03120023

02900040

02900022

02900039

02900104

03210006

00680027

02900020

02900078

04440007

kit of 8 tubes for pump

protection foils

tube-tube
connector


LED-handpiece

mobile
handpiece
holder

handpiece complete
with cord

.....> INSERT SETS

Set of inserts for different applications

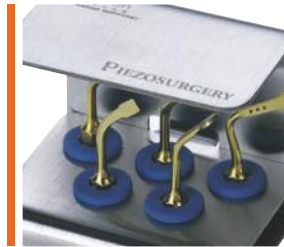

→ BASIC SET



→ OSTEOTOMY SET



→ INSERT SETS

Product	Ref. no.
basic set equipped with inserts OT2, OT7, OP1, OP3 and EX1 description of the single inserts can be found on page 26, 27, 30 and 34	01520001
osteotomy set equipped with inserts OT7, OT7S-4, OT7S-3, OT8R and OT8L description of the single inserts can be found on page 27	01520010

→ INSERT SETS

Product Ref. no.

sinus lift lateral set* 01520023

equipped with inserts
SLC, SLO-H, SLS, SLE1 and SLE2
description of the single inserts can
be found on page 28, 31 and 32

piezo lift set 01520024

equipped with inserts
PL1, PL2 and PL3
description of the single inserts can
be found on page 29

* developed in collaboration with Prof. Tomaso Vercellotti
and Dr. Philippe Russe

→ SINUS LIFT LATERAL SET



→ PIEZO LIFT SET



→ SINUS LIFT SET



→ IMPLANT PREP SET PRO



→ INSERT SETS

Product	Ref. no.
sinus lift set equipped with inserts OT1, OT5, EL1, EL2 and EL3 description of the single inserts can be found on page 26 and 32	01520002
implant prep set pro equipped with inserts IM1S, IM2A, IM2P, IM3A, IM3P, OT4, IM4A, IM4P, P2-3, P3-4, 3 PINS IM1S, 3 PINS 2-2,4 description of the single inserts can be found on page 23, 24, 25 and 26	01520011

→ INSERT SETS

Product	Ref. no.
extraction set equipped with inserts EX1, EX2, EX3, PS2 and PS6 description of the single inserts can be found on page 34 and 35	01520003
third molar extraction set* equipped with inserts EXL1, EXL2, EXL3, EX1, SLO-H description of the single inserts can be found on page 34	01520025
explantation set equipped with inserts EXP3-R, EXP3-L, EXP4-R, EXP4-L description of the single inserts can be found on page 34	01520021

* developed in collaboration with Dr. Fabrizio Fontanella

→ EXTRACTION SET



→ THIRD MOLAR EXTRACTION SET




→ EXPLANTATION SET



→ RESECTIVE PERIO SET



→ RETRO SURGICAL SET

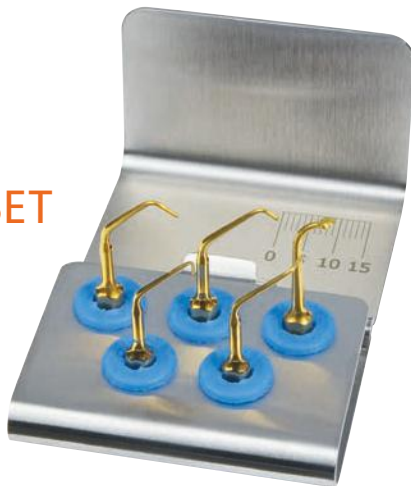

→ INSERT SETS

Product	Ref. no.
resective perio set * equipped with inserts OT13, OT14, OP5A, OP8** and OP9 description of the single inserts can be found on page 28, 30 and 31	01520022
retro surgical set equipped with inserts OP7, EN1, EN3, EN5R, EN5L description of the single inserts can be found on page 31 and 33	01520005

* developed in collaboration with Professor Leonardo Trombelli from the University of Ferrara

** only available for PIEZOSURGERY® 3, PIEZOSURGERY® touch and PIEZOSURGERY white

INDICATIONS

Surgical choice: all

→ SINUS LIFT
CRESTAL APPROACH



→ PIEZO-LIFT

→ PL1 → PL2 → PL3

→ SINUS LIFT
LATERAL APPROACH



→ STANDARD

→ OPTIONAL

→ SLC → SLO-H → SLS → SLE1 → SLE2 → OP3 → OT1 → OT1A → OT5 → OT5A → OT5B → EL1 → EL2 → EL3

→ IMPLANT SITE
PREPARATION



→ STANDARD

→ IM1S → IM2A → IM3A → IM4A → IM2P → IM3P → IM4P → OT4 → P2-3 → P3-4

→ OPTIONAL

→ IM1 AL → IM2A-15 → IM2.8A → IM3A-15 → IM3.4A → IM2P-15 → IM2.8P → IM3P-15 → IM3.4P

→ RIDGE EXPANSION




→ STANDARD

→ OPTIONAL

→ OT7 → OT4 → OP5 → OT2 → OT7A → OT7S-4 → OT7S-3 → OT12 → OT12S → OT7-20

→ PERIOSTEUM
PREPARATION

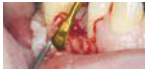

→ STANDARD

→ PR1 → PR2

→ EXTRACTIONS



→ STANDARD

→ EX1 → EX2 → EX3 → PS2

→ THIRD MOLAR
EXTRACTION



→ STANDARD

→ OPTIONAL

→ EXL1 → EXL2 → EXL3 → EX1 → SLO-H

→ STANDARD

→ EXP3-R → EXP3-L → EXP4-R → EXP4-L

→ STANDARD

→ OPTIONAL

→ OT7 → OP5 → OT8R → OT8L → OT6 → OT7A → OT7S-4 → OT7S-3 → OT12 → OT12S → OT7-20

→ STANDARD → OPTIONAL

→ OP3 → OP1 → OP2 → OP3A

→ STANDARD

→ OPTIONAL

→ OP7 → PS2 → EN1 → EN2 → EN3 → EN4 → EN5R → EN5L → EN6R → EN6L → OP3

→ STANDARD

→ OT1 → OT5

→ STANDARD → OPTIONAL

OT7S-4 → OT7S-3 → OT7 → OT7A

→ STANDARD

→ OP5A → OP8* → OP9 → OT13 → OT14 → OP5 → OP3 → PS2 → PP1 → ICP + IC1

→ OPTIONAL

→ PS1 → PS6 → PP10 → PP11 → PP12 → OP2 → OP3A → OP4 → OP6 → OP6A

→ STANDARD

→ DB2 → CROWN PREP TIPS



→ EXPLANTATION



→ BONE BLOCK GRAFTING



→ BONE CHIP GRAFTING/
BONE MODELING



→ ENDODONTICS



→ OSTEOTOMY
CLOSE TO NERVES



→ CORTICOTOMY
TECHNIQUE



→ PERIODONTAL SURGERY
AND IMPLANT CLEANING



→ PROSTHETICS

* only available for PIEZOSURGERY® 3, PIEZOSURGERY® touch and PIEZOSURGERY white

→ TRAYS

Do you like boxing? Don't miss them then.

→ PIEZOSURGERY® INSERT BOX

- fully autoclavable
- perfect storage of all PIEZOSURGERY® inserts
- light, easy to use, compact and roomy

Ref. no. 02900111

→ INSERT TRAY

- stainless steel tray with depth markings
- for storage

Ref. no. 02900037



PIEZOSURGERY® insert box

→ INSERTS

For every clinical needs, the proper insert.

→ IMPLANT PREP



IM1S



PIN IM1S



PIN 2-2.4



PIN 2-3



IM1-AL



PIN IM1-AL



PIN 2-2.4-AL



PIN 2-3-AL

→ IMPLANT PREP INSERTS


Product	Ref. no.
insert IM1S initial pilot osteotomy diamond grain size 30 µm	03510019
PIN IM1S	03740007
PIN 2-2.4	03740003
PIN 2-3 alignment pins, dedicated to IM1S, IM2 (Ø 2), IM3 (Ø 3), OT4 (Ø 2.4), IM1 SP	03740001
insert IM1 AL initial pilot osteotomy diamond grain size 30 µm	03510010
PIN IM1 AL	03740004
PIN 2-2.4 AL	03740005
PIN 2-3 AL 30 mm alignment pins, for extraction sites dedicated to M1 AL, IM2A (Ø 2), OT4 (Ø 2.4), IM3A (Ø 3)	03740006

IMPLANT PREP INSERTS


Product	Ref. no.
insert IM2A	03510001
insert IM2A-15	03510015
insert IM2P	03510002
insert IM2P-15	03510016
pilot osteotomy in anterior and posterior region	
insert IM2.8A	03510025
insert IM2.8P	03510026
insert IM3A	03510003
insert IM3A-15	03510017

to enlarge or to finalize the implant site preparation;
insert with double irrigation to avoid overheating

IMPLANT PREP

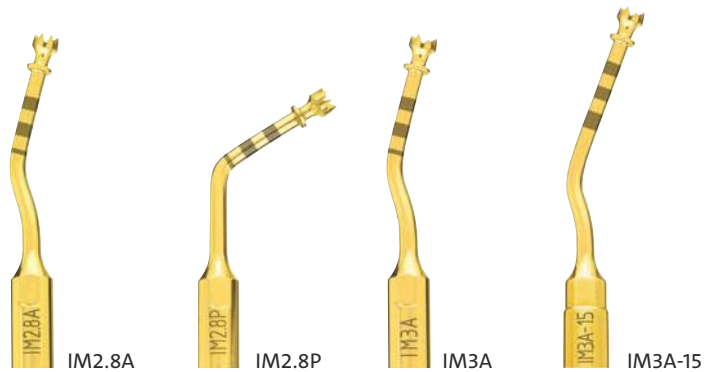

→ IMPLANT PREP

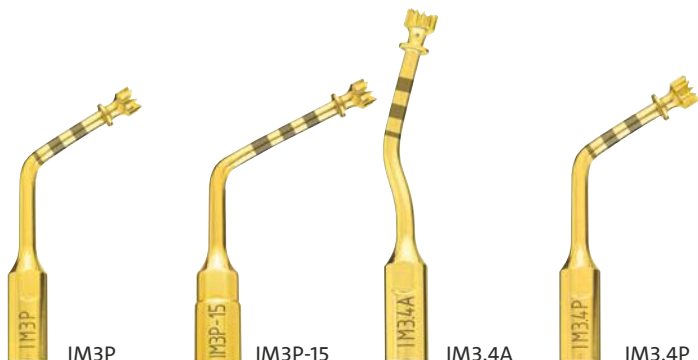

IM3P

IM3P-15

IM3.4A

IM3.4P



IM4A

IM4P

P2-3

P3-4

→ IMPLANT PREP INSERTS

Product	Ref. no.
insert IM3P	03510004
insert IM3P-15	03510018
insert IM3.4A	03510027
insert IM3.4P	03510028
to enlarge or to finalize the implant site preparation; insert with double irrigation to avoid overheating	
insert IM4A	03510005
insert IM4P	03510006
to finalize the implant site preparation; insert with double irrigation to avoid overheating	
insert P2-3	03510020
insert P3-4	03510021
to optimize concentricity of implant site preparation between $\varnothing 2$ and $\varnothing 3$ mm diamond grain size 150 μm	

OSTEOTOMY INSERTS

Product	Ref. no.
insert OT1	03370001
insert OT1A	03370008

to finalize the osteotomy in proximity of soft tissues
(for example: sinus membrane, vessel, alveolar nerve)
OT1 diamond grain size 91 µm
OT1A diamond grain size 126 µm

insert OT2 03370002
osteotome of great precision in anatomically thin
structures (for example, ridge expansion, interdental
corticotomies, non traumatic nasal spine)

insert OT4 03370004
to correct pilot osteotomy axis, to finalize the implant
site preparation close to the alveolar nerve; sinus crestal
approach technique
diamond grain size 151 µm

insert OT5 03370005
insert OT5A 03370009
insert OT5B 03370010

non-traumatic, to finalize the osteotomy or osteoplasty
on thin bone and/or near delicate anatomic structures
OT5 Ø 1,7 - diamond grain size 91 µm
OT5A Ø 1,9 - diamond grain size 91 µm
OT5B Ø 1,7 - diamond grain size 126 µm

OSTEOTOMY



OSTEOTOMY



OT6
5 teeth,
thickness 0.75 mm

OT7A
5 teeth,
thickness 0.75 mm

OT7
5 teeth,
thickness 0.55 mm

OT7-20
5 teeth,
thickness 0.6 mm



OT7S-4
4 teeth,
thickness 0.35 mm



OT7S-3
3 teeth,
thickness 0.35 mm



OT8R
thickness 0.6 mm



OT8L
thickness 0.6 mm

OSTEOTOMY INSERTS

Product	Ref. no.
insert OT6	03370006
osteotomy of large bone sections during maxillofacial surgery	
insert OT7	03370007
insert OT7A	03370011
all the osteotomy technique in maxilla and mandible	
– ridge expansion	
– corticotomy technique	
– bone block grafting	
insert OT7-20*	03370017
ridge expansion technique, bone block grafting (from chin/mandible ramus), Le Fort I osteotomy, bilateral sagittal split osteotomy	
insert OT7S-4	03370014
insert OT7S-3	03370015
very thin osteotomy, corticotomy for orthodontic microsurgery technique, root separation in dental extraction technique and periodontal surgery	
insert OT8R	03370012
insert OT8L	03370013
all the osteotomy technique in maxilla and mandible	
– bone block grafting	

* tighten up this insert only with torque wrench K8

OSTEOTOMY INSERTS

Product Ref. no.

insert OT9 03370016

to finalize the osteotomy or osteoplasty near delicate anatomic structures; removal of the sinus basal cortex; insert with double irrigation to avoid overheating

insert OT11 03370019

removal of the sinus basal cortex and elevation of the sinus membrane during crestal sinus lift procedures

insert OT12 03370020

insert OT12S 03370021

dedicated to the osteotomy techniques where the standard piezoelectric bone saws may have a difficult surgical approach

insert OT13 03370022

insert OT14 03370023

non-traumatic bone defect preparation on buccal and/or palatal cortical plates in resective periodontal surgery


OT13 \varnothing 1.7 - diamond grain size 150 μ m

OT14 \varnothing 2.3 - diamond grain size 150 μ m

insert SLO-H 03370025

bone window antrostomy in sinus lift procedure by vestibular approach

OSTEOTOMY

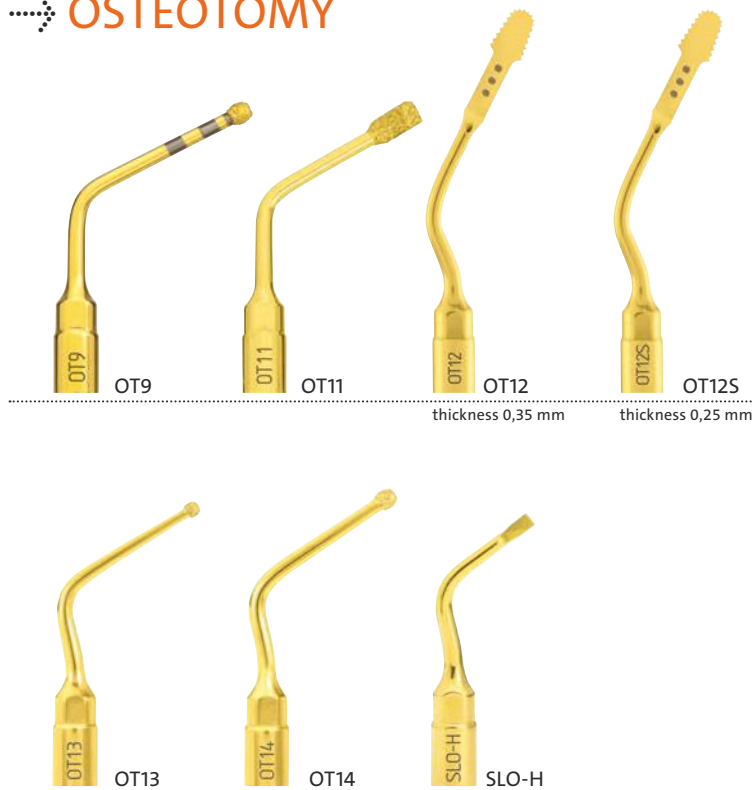

→ OSTEOTOMY



→ OSTEOTOMY INSERTS

Product	Ref. no.
insert PL1	03370026
sinus floor reaching and safe bony ring removal during crestal sinus lift procedure (PIEZO-LIFT)	
insert PL2	03370027
sinus floor consumption and initial membrane elevation during crestal sinus lift procedure (PIEZO-LIFT)	
insert PL3	03370028
removal of the sinus basal cortex and elevation of the sinus membrane during crestal sinus lift procedure (PIEZO-LIFT)	

OSTEOPLASTY INSERTS

Product	Ref. no.
insert OP1 bone remodelling and harvesting of bone chips	03380001
insert OP2 bone osteotomy, remodelling and removal of inflammatory tissue	03380002
insert OP3 insert OP3A periodontal osteotomy, crown lengthening, bone chips harvesting, inflammatory tissue removal (cyst, etc.)	03380003 03380011
insert OP4 interproximal osteoplasty and root planing diamond grain size 126 µm	03380004
insert OP5 insert OP5A root debridement and root planing during resective and regenerative periodontal surgery OP5 diamond grain size 30 µm OP5A diamond grain size 90 µm	03380005 03380012
insert OP6 root preparation in periodontal surgery OP6 diamond grain size 126 µm	03380006

OSTEOPLASTY



→ OSTEOPLASTY



→ OSTEOPLASTY INSERTS

Product	Ref. no.
insert OP7	03380007
peri-apical maxillary bone osteotomy access, inflammatory tissue removal	
insert OP8*	03380008
insert OP9	03380009
Interproximal osteoplasty without damaging adjacent root surfaces	
OP8* cutting edge from 0.7 - 1.3 mm	
OP9 cutting edge from 1 - 2 mm	
insert SLC	03380013
bone consuming during lateral sinus lift procedures	

* only available for PIEZOSURGERY® 3, PIEZOSURGERY® touch
and PIEZOSURGERY white

→ ELEVATION INSERTS

Product	Ref. no.
insert EL1 separation of the sinus membrane, 2 mm around the frame of bony window	03390001
insert EL2	03390002
insert EL3 separation of the sinus membrane in internal zones	03390003
insert SLS separation of the sinus membrane during lateral sinus lift procedures	03390006
insert SLE1 elevation of the sinus membrane from the sinus floor	03390007
insert SLE2 elevation of the sinus membrane from the palatal sinus wall	03390008

→ PERIOSTEUM INSERTS

insert PR1	03390004
insert PR2 subperiosteal preparation	03390005

→ ELEVATION



→ PERIOSTEUM



ENDO RETRO



EN1



EN2



EN3



EN4



EN5R



EN5L



EN6R



EN6L

ENDO RETRO INSERTS

Product	Ref. no.
insert EN1 diamond-coated endo apical debrider 3 mm, apical root debridement diamond grain size 30 µm	02170001
insert EN2 smooth endo apical debrider 3 mm, gentle apical root debridement	02170002
insert EN3 diamond-coated endo apical debrider 2.2 mm, apical root debridement diamond grain size 30 µm	02170005
insert EN4 smooth endo apical debrider 2.2 mm, gentle apical root debridement	02170006
insert EN5R	02170007
insert EN5L right and left angled, diamond-coated endo apical debrider 2.2 mm, apical root debridement diamond grain size 30 µm	02170008
insert EN6R	02170009
insert EN6L right and left angled, smooth endo apical debrider 2.2 mm, apical root debridement	02170010

→ EXTRACTION INSERTS

Product	Ref. no.
insert EX1 to cut off the ankylosis, root fraction techniques	03400001
insert EX2 analogue to EX1 in the posterior regions	03400002
insert EX3 analogue to EX1 in the posterior regions	03400003
insert EXL1 for luxation and simultaneous root extraction in a single manoeuvre	03400008
insert EXL2 analogue to EXL1. It has a smaller radius but makes it possible to apply greater force	03400009
insert EXL3 for alveolar debridement and/or removing radicular fragments from the extractive alveolus	03400010

→ EXPLANTATION INSERTS

insert EXP3-R	03400004
insert EXP3-L	03400005
insert EXP4-R	03400006
insert EXP4-L	03400007

explantations

→ EXTRACTION



→ EXPLANTATION





PS1



PS2



PS6



PP1



PP10



PP11



PP12

PERIODONTAL INSERTS

Product	Ref. no.
insert PS1 root scaling	03180001
insert PS2 scaling and inflammatory tissue removal fractured root apex extraction	03180003
insert PS6 root scaling	03180005
insert PP1 root planing	03180002
insert PP10 gentle subgingival concrements removal	03180006
insert PP11 gentle subgingival concrements removal	03180007
insert PP12 gentle subgingival concrements removal	03180008

→ IMPLANT CLEANING



→ IMPLANT CLEANING INSERT

Product	Ref. no.
Implant cleaning set P (consists of a tipholder ICP and 5 tips IC1)	01520020
tipholder ICP	03570003
tip IC1 (5 pieces)	02900112

efficient and soft cleaning of implants and restorations; implant cleaning tip IC1 can be easily tightened without any special tools and is fully autoclavable


→ CROWN PREP



key AB1



dynamometric wrench K7

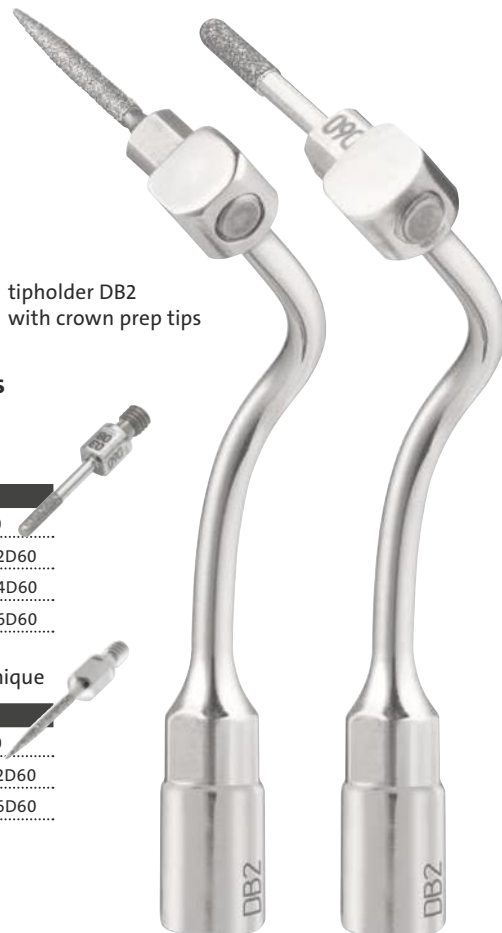

tipholder DB2 with crown prep tips

→ Available lengths and diameters of crown prep tips

Preparation with chamfer technique

→ Length 10 mm → Diamond coating

diameter	D120	D90	D60
1.2 mm		TA12D90	TA12D60
1.4 mm	TA14D120	TA14D90	TA14D60
1.6 mm	TA16D120	TA16D90	TA16D60

Preparation with vertical preparation technique

→ Length 10 mm → Diamond coating

diameter	D90	D60
1.2 mm	TF12D90	TF12D60
1.6 mm	TF16D90	TF16D60

→ CROWN PREP TIPHOLDER FOR PIEZOSURGERY® HANDPIECES

Product	Ref. no.
crown prep set PS equipped with DB2, AB1 and K7 tipholder DB2	02160010 03570002

→ CROWN PREP TIPHOLDER ACCESSORIES

key AB1	03580001
dynamometric wrench K7	02900081

→ CROWN PREP TIPS FOR CHAMFER T.

crown prep tip TA12D60	03590001
crown prep tip TA14D60	03590002
crown prep tip TA16D60	03590003
crown prep tip TA12D90	03590004
crown prep tip TA14D90	03590005
crown prep tip TA16D90	03590006
crown prep tip TA14D120	03590007
crown prep tip TA16D120	03590008

→ CROWN PREP TIPS FOR VERTICAL T.


crown prep tip TF12D60	03590010
crown prep tip TF16D60	03590011
crown prep tip TF12D90	03590012
crown prep tip TF16D90	03590013



→ MULTIPIEZO

Looking beautiful is looking really beautiful.


RESTORATIVE MODE
PULSE MODE

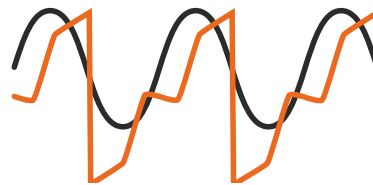

→ Pulse Mode

Pulse Mode (level 6 of restorative function) transforms sine wave-shaped ultrasound oscillations to a new profile, characterized by specific power peaks, which results in an optimized performance for extractions and prosthetics.

STANDARD MODE
SOFT MODE

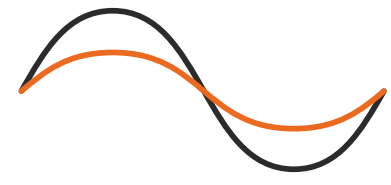

→ SOFT MODE

– special modulation for gentle but still powerful insert movement
– maximum comfort for sensitive patients

→ MULTIPIEZO

Equipment set

- 1 MULTIPIEZO
- 1 LED ultrasonic handpiece slim
- 2 inserts: scaling insert S1-S and perio insert P3
- 2 bottles, 500 ml
- 2 safety bottle caps

Ref. no.

05060023



→ TECHNICAL DATA

sizes (l x w x h in mm):

320 x 230 x 145

weight:

2.4 kg

operating frequency:

24.000 – 36.000 Hz

power types:

3 functions, each with seven power levels, SOFT MODE

liquid supply:

via removable bottle

liquid quantity:

adjustable in 7 levels

liquid flow rate:

between 0 and 28 ml

warranty:

2 years, handpiece 1 year

→ MICROPIEZO

A new era in ultrasound devices.

The advanced entry level ultrasound scaler with three functions:

- scaling and perio treatments
- endodontic treatments
- prosthetic treatments

Patients benefit from gentle ultrasound scaling in SOFT MODE.



→ MICROPIEZO

- 1 MICROPIEZO
- 1 LED ultrasonic handpiece slim
- 2 inserts: scaling insert S1-S and perio insert P3
- 1 liquid container, 500 ml
- 1 footpedal

Ref. no.

05080013





→ TECHNICAL DATA

sizes (l x w x h in mm):

161 x 214 x 171

weight:

0.94 kg

operating frequency:

24.000 – 36.000 Hz

power types:

3 functions (endo, perio-scaling, prosthetics),
SOFT MODE, each with seven power levels

liquid supply:

via removable bottle

liquid quantity:

adjustable in 7 levels

liquid flow rate:

between 0 and 33 ml

→ TECHNICAL DATA

operating frequency:

24.000 – 36.000 Hz

power types:

adjustable via the dental chair

power supply:

24 VAC, 50-60 Hz / 32 VDC

water supply:

via the dental unit

warranty:

2 years, handpiece 1 year

→ COMPACT PIEZO LED

Equipment set

1 electronic module

1 LED ultrasonic handpiece slim

2 inserts: scaling insert S1-S and perio insert P3

Ref. no. 05090008

→ COMPACT PIEZO P2K

Equipment set

1 electronic module

1 ultrasonic handpiece

2 inserts: scaling insert S1-S and perio insert P3

Ref. no. 05090003

COMPACT PIEZO LED COMPACT PIEZO P2K

Compact efficiency or directly built-in.



Including feedback-system
for constant power.

SPARE PARTS

High quality spare parts for your ultrasound device.



DSC adaptor for ultrasound inserts

scaler insert adaptor

ADAPTORS FOR THERMODISINFECTION

Product

- DSC Adaptor (for scaler handpiece irrigation line) *
- spare filter for DSC Adaptor *
- scaler insert Adaptor *

Ref. no.


- 04610007
- 04590006
- 04610009

COMMON SPARE PARTS

- dynamometric torque wrench K10
- dynamometric torque wrench K6
- dynamometric wrench K7 (for tipholder DB1)
- water filter
- insert tray
- universal cart

- 02900137-001
- 02900074
- 02900081
- 00420004
- 02900072
- 03540011

* not disposable

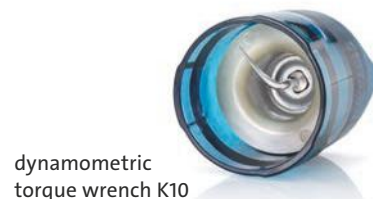

dynamometric torque wrench K10



dynamometric torque wrenches K6 and K7



water filter



safety
bottle cap
(blue)



liquid container (MICROPIEZO),
incl. cap, 500ml



plastic bottle
(MULTIPIEZO)
with complete
bottle cap

→ SPARE PARTS ULTRASOUND UNITS WITH LED HANDPIECE

Product

Ref. no.

LED ultrasonic handpiece slim

03120238

kit light guide, metal cone and blue O-ring

for LED ultrasonic handpiece slim

02900146

LED ultrasonic handpiece

03120142

handpiece LED-cone

03020158

handpiece metal cone

03020171

→ SPARE PARTS ULTRASOUND UNITS WITHOUT LED HANDPIECE

ultrasonic handpiece

03120016

handpiece metal cone

00710034

→ BOTTLES AND BOTTLE CAPS FOR PIEZO SMART AND MULTIPIEZO RANGE

plastic bottle (piezo smart), 300 ml,

with blue complete bottle cap

03020066

plastic bottle (piezo smart), 300 ml,

with grey complete bottle cap

03020067

plastic bottle (piezo smart), 300 ml, without cap

03230007

plastic bottle (MULTIPIEZO), 500 ml, without cap

03150102

complete bottle cap, blue

03020081

complete bottle cap, grey

03020082

safety bottle cap, grey

01950015

safety bottle cap, blue

01950014

→ SPARE PARTS MULTIEZEO RANGE

Product	Ref. no.
O-ring set for bottle holder (6 pieces)	03020175
peristaltic pump (white)	03210001
peristaltic pump (black - from serial no. 106005993)	03210009
footswitch	02900114

→ SPARE PARTS MICROPIEZO

LED ultrasonic handpiece slim	03120238
kit light guide, metal cone and blue O-ring for LED ultrasonic handpiece slim	02900146
liquid container (incl. cap, 500ml)	02900145
footpedal	02900144
pump	03210009
power adapter	03160024
kit O-ring + spring for liquid container	02900147

→ SPARE PARTS PIEZO SMART

O-ring set for bottle holder (3 pieces)	03020068
peristaltic pump	03210007
footswitch	02900055
power adapter piezo smart	03020105

→ SPARE PARTS MICROPIEZO S

footswitch	02900001
water filter	00420004



LED ultrasonic handpiece slim



LED ultrasonic handpiece



ultrasonic handpiece


→ INSERT SETS

An insert set for every indication.



→ K10 – protects your handpiece and inserts

- sterilisable dynamometric torque wrench
- supplied with inserts

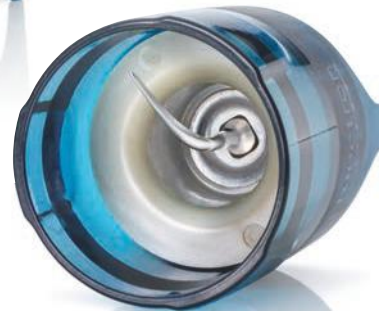

→ SET MAINTENANCE



→ SET SCALING



→ INSERT SETS

Product

Ref. no.

set maintenance

02160035-001

equipped with inserts
S1, S2, S5, P2, P3 and
5 torque wrenches K10

description of the single inserts
can be found on page 50 and 51

set scaling

02160001-001

S1, S2, S3, S4, S5 and
5 torque wrenches K10

description of the single inserts
can be found on page 50

INSERT SETS

Product

Ref. no.

set perio maintenance

02160006-001

equipped with inserts
P10, P11, P12, P13, P14 and
5 torque wrenches K10
description of the single inserts
can be found on page 52

set perio treatment

02160036-001

equipped with inserts
P2, P3, P15, P16R, P16L and
5 torque wrenches K10
description of the single inserts
can be found on page 51 and 52

SET PERIO MAINTENANCE



SET PERIO TREATMENT



→ SET IMPLANT CLEANING



→ SET ENDO REVISION



→ INSERT SETS

Product

Ref. no.

implant cleaning set S

02160033-001

equipped with inserts

1 insert ICS

1 torque wrench K10

5 Implant cleaning tips IC1

description of the single inserts
can be found on page 53

set endo revision

02160008-001

equipped with inserts

ER1, ER2, ER3, ER4, ER5 and

5 torque wrenches K10

description of the single inserts
can be found on page 55



→ every insert
is supplied with its
own torque wrench


→ SCALING INSERTS

Product	Ref. no.
insert S1 universal curette with semicircular diameter for considerable tartar removal	02960001-001
insert S1-S slim universal curette with triangled slightly curved surface for removal of supragingival tartar, subgingival concretions and biofilm	02960009-001
insert S2 universal curette with triangled, slightly curved surface, efficient in the interdental spaces and posterior surfaces	02960002-001
insert S3 flat, with rounded edges for considerable supragingival tartar removal	02960003-001
insert S4 universal curette with 45° angled, triangled and slightly curved surface for tartar removal on mesial and distal interdental surfaces in the posterior area	02960004-001
insert S5 similar shape to S1 insert but longer and thinner for gentle supra- and subgingival tartar removal and for gingivitis	02960005-001
insert S6 contra-angled universal curette with semicircular diameter, powerful insert (twice the power of S1 insert) for considerable tartar removal	02960006-001

→ INSERTS

For a wide range of applications

→ SCALING

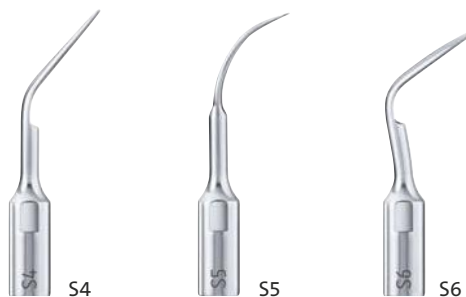

→ PERIO UNIVERSAL



→ every insert
is supplied with
its own torque wrench

→ PERIO UNIVERSAL INSERTS

Product	Ref. no.
insert P1 long, straight working tip with a circular surface	02970001-001
insert P2 curved working tip with a circular surface	02970002-001
insert P3 slim insert smooth like a probe	02970003-001
insert P4 extra long, straight working tip with a circular surface	02970004-001



→ every insert
is supplied with its
own torque wrench

→ PERIO ANATOMIC INSERTS

Product	Ref. no.
insert P10 extra long, curved working tip with a circular surface	03080001-001
insert P11 right (inclination of 15°) curved insert with round tip	03080002-001
insert P12 left (inclination of 15°) curved insert with round tip	03080003-001
insert P13 right (inclination of 15°) curved insert with spherical (Ø = 0.8 mm) round tip	03080004-001
insert P14 left (inclination of 15°) curved insert with spherical (Ø = 0.8 mm) round tip	03080005-001
insert P15 universal curette shape with round surface	03080006-001
insert P16R right angled periodontal curette	03080007-001
insert P16L left angled periodontal curette	03080008-001

→ PERIO ANATOMIC



→ IMPLANT CLEANING



→ implant cleaning tip IC1 can be easily tightened without any special tools and is fully autoclavable



→ every insert is supplied with its own torque wrench

→ IMPLANT CLEANING INSERT

Product	Ref. no.
Implant cleaning set S (consists of a tipholder ICS and 5 tips IC1)	02160033-001
tipholder ICS	03570004-001
tip IC1 (5 pieces)	02900112
efficient and soft cleaning of implants and restorations	



→ every insert
is supplied with its
own torque wrench

→ ENDO INSERTS

Product	Ref. no.
insert E1 120° file holder for endo files and inserts with a shaft of 0,8 mm	02350001-001
insert E2 90° file holder for endo files and inserts with a shaft of 0,8 mm	02350002-001

→ SPARE PARTS

key K1 for file holders E1 and E2	02270003
key for screw of file holders E1 and E2	02270004
screw for file holder E1	01260005
screw for file holder E2	01260003

→ ENDO




E1



E2

→ ENDO REVISION

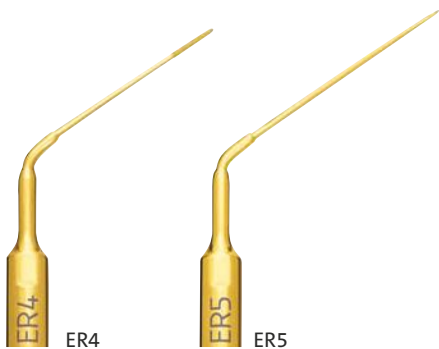

→ every insert is supplied with its own torque wrench

→ ENDO REVISION INSERTS

Product	Ref. no.
insert ER1 angled shaft with small, spherical (\varnothing 1.7 mm) and D30 diamond coated tip	03450001-001
insert ER2 angled insert with 0.7 mm \varnothing , working length up to 10 mm, the last 5 mm D30 diamond coated	03450002-001
insert ER3 angled insert with 0.6 mm \varnothing , working length up to 10 mm, the last 5 mm D30 diamond coated	03450003-001
insert ER4 angled insert with 0.6 mm \varnothing , working length up to 20 mm, the last 5 mm D30 diamond coated	03450004-001
insert ER5 angled insert with 0.5 mm \varnothing , working length up to 24 mm	03450005-001



→ every insert
is supplied with its
own torque wrench

→ ENDO RETRO INSERTS

Product Ref. no.
insert RD3 03050006-001
angled shaft, with 90° angled working tip
and fine diamond coating (D30)

insert RD4 03050007-001
angled shaft, with 90° right angled working tip
and fine diamond coating (D30)

insert RD5 03050008-001
angled shaft, with 90° left angled working tip
and fine diamond coating (D30)




→ equipped with
K6 torque wrench

→ EXTRACTION INSERTS

Product Ref. no.
insert ME1* 03130001-001
insert ME2* 03130002-001
insert ME3* 03130003-001

The specially shaped inserts ME1, ME2 and ME3 – with a blade thickness of only 0.35 mm – perfectly fulfill the anatomical demands of non surgical extractions. They offer maximum intra-operative sensitivity, allowing a great preservation of alveolar bone and periodontal tissues. Furthermore, the physiological solution guarantees maximum visibility thanks to the cavitation effect which leads to a temporary hemostasis of blood vessels.

→ ENDO RETRO

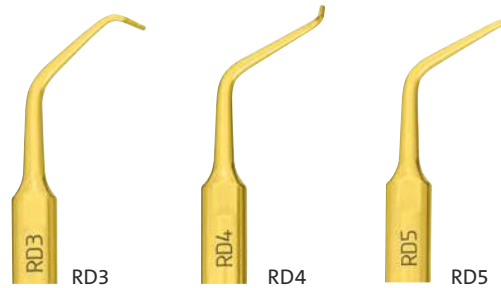

→ EXTRACTION*




→ CROWN PREP



key AB1



dynamometric wrench K7

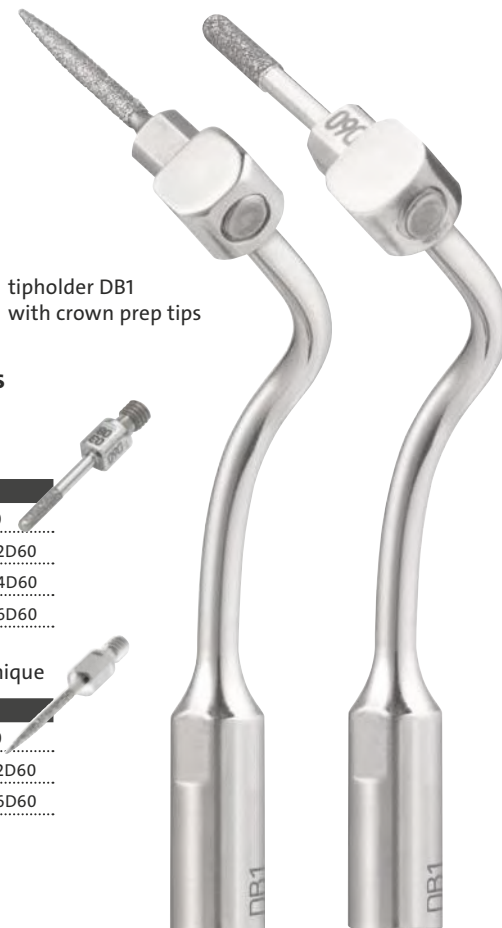

tipholder DB1 with crown prep tips

→ Available lengths and diameters of crown prep tips

Preparation with chamfer technique

→ Length 10 mm → Diamond coating

diameter	D120	D90	D60
1.2 mm		TA12D90	TA12D60
1.4 mm	TA14D120	TA14D90	TA14D60
1.6 mm	TA16D120	TA16D90	TA16D60

Preparation with vertical preparation technique

→ Length 10 mm → Diamond coating

diameter	D90	D60
1.2 mm	TF12D90	TF12D60
1.6 mm	TF16D90	TF16D60

→ CROWN PREP TIPHOLDER FOR ULTRASOUND SCALER HANDPIECES

Product	Ref. no.
crown prep set scaler equipped with DB1, AB1 and K7 tipholder DB1	02160009 03570001

→ CROWN PREP TIPHOLDER ACCESSORIES

key AB1	03580001
dynamometric wrench K7	02900081

→ CROWN PREP TIPS FOR CHAMFER T.

crown prep tip TA12D60	03590001
crown prep tip TA14D60	03590002
crown prep tip TA16D60	03590003
crown prep tip TA12D90	03590004
crown prep tip TA14D90	03590005
crown prep tip TA16D90	03590006
crown prep tip TA14D120	03590007
crown prep tip TA16D120	03590008

→ CROWN PREP TIPS FOR VERTICAL T.

crown prep tip TF12D60	03590010
crown prep tip TF16D60	03590011
crown prep tip TF12D90	03590012
crown prep tip TF16D90	03590013

COMBI touch

Ultrasound and air-polishing in one unit.

- flexible supra- and subgingival air-polishing
- glass touch screen, easy to clean
- removable powder containers
- SOFT MODE for ultra-gentle scaling




COMBI touch PERIO

Equipment set

- 1 LED ultrasonic handpiece slim
- 2 inserts: scaling insert S1-S and perio insert P3
- 1 bottle (500ml)
- 1 grey safety bottle cap
- 1 air-polishing handpiece 120°
- 1 air-polishing handpiece PERIO
- 1 pack Subgingival Perio Tips (40 pcs)
- 1 wrench k9
- 2 bottles prophylaxis powder sensitive+ - glycine powder (160g)

Ref. no.

05020015

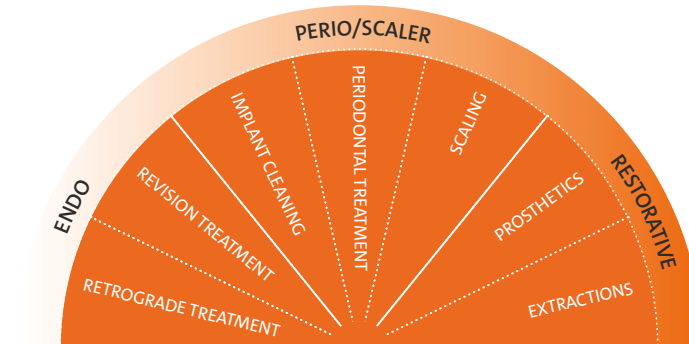

- wide range of ultrasound indications



→ TECHNICAL DATA

sizes (l x w x h in mm): 410 x 260 x 145

weight: 4.8 kg

water supply:

working pressure from 1 to 6 bars

air supply:

input pressure between 4 and 8 bars

ultrasound part:

operating frequency:

24 KHz to 36 KHz

power mode:

"endo", "perio/scaler", "restorative", "SOFT MODE"

power levels from 1 to 5

irrigation:

tap water line – stepless adjustment

peristaltic pump – adjustable on the

touch screen: 7 flow levels from 1

(approx. 5 ml/min) to 7 (approx. 30 ml/min)

possibility to switch off the irrigation in the

power modes "endo", "perio/scaler" and

"restorative" at power levels from 1 to 4

polishing part:

polisher function:

selectable via touch screen: "Prophy" or "Perio"

irrigation:

continuous regulation via knob,

water heated via heating element.

warranty:

2 years, handpiece 1 year

7 different turbine adapters

- Bien-Air UNIFIX® and UNIFIX® L
- KaVo MULTIFLEX® and MULTIFLEX® lux
- SIRONA® R and B
- W&H® ROTO QUICK® and ROTO QUICK® lux
- NSK® MACHLITE® and PHATELUS®
- Midwest
- Borden

- UNIFIX® is a registered trademark of Bien-Air (BIEN-AIR®) Holding SA
- MULTIFLEX® is a registered trademark of Kaltenbach & Voigt (KaVo®) GmbH
- SIRONA® is a registered trademark of SIRONA DENTAL SYSTEMS GMBH
- ROTO QUICK® is a registered trademark of W&H® Deutschland GmbH
- MACHLITE® and PHATELUS® are registered trademarks of Nakanishi Corporation and Nakanishi Shika Kikai (NSK®) Ltd

STARJET BASIC

Equipment set

1 STARJET

1 spray head 120°

2 bottles prophylaxis powder sensitive+ - glycine powder (160g)

Product

STARJET basic (BienAir)
STARJET basic (KaVo)
STARJET basic (Sirona)
STARJET basic (W&H)
STARJET basic (Midwest)
STARJET basic (Borden)
STARJET basic (NSK)

Ref. no.

05150074
05150075
05150076
05150077
05150078
05150079
05150080

STARJET

The 2 face air-polisher.



One device, two indications – easy switch from supra- to subgingival air-polishing with a PROPHY to PERIO turn.

- maximum flexibility as it can be used with powders with different granulometries
- optimal access thanks to 3 available spray heads – 120°, PERIO and 90° optional
- permanent control of the powder level thanks to the transparent powder chamber
- easy powder filling thanks to lateral opening



→ TECHNICAL DATA

sizes (l x w x h in mm):

200 x 90 x 50

water supply:

0 – 3.0 bar

air supply:

2.0 – 3.0 bar

warranty:

1 year

→ STARJET

Equipment set

1 STARJET

1 spray head 120°

1 spray head PERIO;

Subgingival Perio Tips (40 pcs)

1 wrench K9

2 bottles prophylaxis powder sensitive+ - glycine powder (160g)

Product

Ref. no.

STARJET (BienAir)

05150074

STARJET (KaVo)

05150075

STARJET (Sirona)

05150076

STARJET (W&H)

05150077

STARJET (Midwest)

05150078

STARJET (Borden)

05150079

STARJET (NSK)

05150080

MECTRON PROPHYLAXIS POWDERS

Higher performance for bright results.

PROPHYLAXIS POWDER

Product	Ref. no.
Rhapsody Flow Lemon* pack with four bottles (300 g each)	03140039
Rhapsody Flow Neutral* pack with four bottles (300 g each)	03140040
Prophylaxis powder sensitive+ pack with two bottles (160 g each)	03140021

* Rhapsody Flow powders are registered and manufactured by Queen Dental GmbH, Neue Straße 67, 99846 Seebach, Germany

Only available in EU countries - Distributor: Mectron Spa

→ RHAPSODY FLOW LEMON

gently and
effectively
removes
supragingival
biofilm and
stains on the
enamel surface



- Sodium bicarbonate-based
- Average particle size 40 µm
- Lemon taste

→ RHAPSODY FLOW NEUTRAL

gently and
effectively
removes
supragingival
biofilm and
stains on the
enamel surface



- Sodium bicarbonate-based
- Average particle size 40 µm
- Neutral taste

WE ♥ PRO
PHYLAXIS

WE ♥ PRO
PHYLAXIS

WE ♥ PRO
PHYLAXIS

→ PROPHYLAXIS
POWDER SENSITIVE⁺

gently removes
supra- and sub-
gingival biofilm



- Glycine-based
- Average particle size 25 µm
- Sweet taste

INDICATIONS	SODIUM BICARBONATE	GLYCINE
Supragingival biofilm removal	●	●
Removal of stains and pigmentations	●	
Subgingival biofilm removal		●
Periodontal pocket detoxification		●
Implant maintenance		●
Peri-implantitis treatment		●
Fixed orthodontic appliances	●	●
Cleaning of prosthesis and restoration materials		●
Restoration maintenance	●	●

→ SIZES:

total (l x w x h in mm):

480 x 460 x 850 mm

upper tray:

480 x 370 in 770 mm height

middle tray:

420 x 310 in 440 mm height

lower tray:

480 x 370 in 65 mm height

weight:

29 kg

→ UNIVERSAL CART

Let's move the devices!

The **universal cart** – maximum mobility for all your table-top devices. All supply lines are hidden in the easily accessible cable room. The power supply is brought at the tray's height, for an easy connection to your devices.



→ UNIVERSAL CART

Product

universal cart

Ref. no.

03540011



water filter



plastic bottle (scaler) –
COMBI touch



LED ultrasonic handpiece slim



LED ultrasonic handpiece



dynamometric
torque wrench K6



dynamometric
torque wrench K10



SPARE PARTS

COMMON SPARE PARTS COMBI touch

Product

Ref. no.

dynamometric torque wrench K10	02900137-001
dynamometric torque wrench K6	02900074
dynamometric torque wrench K7	02900081
plastic bottle (scaler), 500 ml, without cap	03150102
complete bottle cap, blue	03020081
complete bottle cap, grey	03020082
safety bottle cap, grey	01950014
safety bottle cap, blue	01950015
Subgingival Perio Tip (40 pieces)	02900120
wrench K9	02900122
cleaning needles kit 0,4 x 57 mm and 0,8 x 60 mm	02900151
O-ring set COMBI touch	03020206
powder chamber PROPHY	03020194
powder chamber PERIO	03020200
powder chamber cap	03020195
water filter	00420004
water supply hose	02900119
air supply hose	02900012
suitcase COMBI touch	04440019
footswitch	02900114

SPARE PARTS COMBI *touch* WITH SLIM HANDPIECES

Product	Ref. no.
LED ultrasonic handpiece slim	03120238
kit light guide metal cone and blue O-ring for LED ultrasonic handpiece slim	02900146
air-polishing handpiece 120°	03120240
air-polishing handpiece PERIO	03120241
air-polishing handpiece 90°	03120239
cleaning needles kit 0,4 x 57 mm and 0,8 x 60 mm	02900151

SPARE PARTS COMBI *touch* WITH STANDARD HANDPIECES up to serial no. 205010618

LED ultrasonic handpiece	03120224
handpiece LED-cone	03020205
handpiece metal cone	03020171
90° spray head	03340007
120° spray head	03340008
PERIO spray head	03340009
O-ring set COMBI <i>touch</i>	03020206
cleaning needles kit 0,4 x 57 mm and 0,8 x 57 mm	02900149

COMMON SPARE PARTS COMBI AND TURBODENT

dynamometric torque wrench K10	02900137-001
torque wrench K6	02900074
cleaning needle	02900003
O-ring for powder chamber cap	00410018
footswitch	02900001
water filter	00420004
O-ring water coupling	00410024



air-polishing handpiece 120°



air-polishing handpiece PERIO



air-polishing handpiece 90°



90° spray head



120° spray head



PERIO spray head



ultrasonic handpiece



spray head set turbodent p/s



spray head, to easyjet serial no. 080XXXXY



spray head, from easyjet serial no. 080XXXXY



O-ring Set

→ ADAPTORS FOR THERMODISINFECTION

Product	Ref. no.
DSC Adaptor (for scaler handpiece irrigation line) *	04610007
spare filter for DSC Adaptor *	04590006
scaler insert Adaptor *	04610009

→ SPARE PARTS COMBI, up to serial no. 046XXXXY / TURBODENT, up to serial no. 047XXXXY

ultrasonic handpiece	03120016
ultrasonic handpiece cord (grey)	02830022
spray head set turbodent p/s	03340003
air polishing handpiece cord (grey)	03120009
O-ring spray head, small	00410013
O-ring spray head, large	00410014

→ SPARE PARTS COMBI, from serial no. 203XXXXY / TURBODENT, from serial no. 202XXXXY

ultrasonic handpiece cord (white)	02830042
spray head	03340006
spray head set (spray head, silicon protector, long cleaning needle and short cleaning needle)	03020121
air polishing handpiece cord	03120133
silicon protector	03150033

→ SPARE PARTS STARJET

90° spray head	03340015
120° spray head	03340014
PERIO spray head	03340013
Subgingival Perio Tip (40 pieces)	02900120
wrench K9	02900122

* not disposable

→ SPARE PARTS STARJET

Product	Ref. no.
O-ring set STARJET (2 black + blue)	02900134-001
set powder chamber cap + 1 O-ring blue	02900133-001
cleaning kit STARJET	02900136-001
(cleaning needle spray head, cleaning needle STARJET body, cleaning brush)	

→ COMMON SPARE PARTS EASY JET PRO

cleaning needle (long)	02900050
silicon protector for spray head	03150033
O-ring Borden turbine adapter	01970005
O-ring Midwest turbine adapter	01970013

→ SPARE PARTS EASYJET PRO, up to serial no. 080XXYYYY

spray head set easyjet pro	03020061
(spray head, silicon protector, cleaning needle, small and large O-ring)	
O-ring set easyjet pro	02900051
(1 x small and 1 x large O-ring for the spray head, 1 x O-ring for the powder chamber cap)	
powder chamber cap, transparent	03020059
powder chamber cap, blue	03020060

→ SPARE PARTS EASYJET PRO (from serial no. 201XXYYYY) AND EASYJET PERIO

spray head	03340006
spray head set	03020121
(spray head, silicon protector, long cleaning needle and short cleaning needle)	
powder chamber cap, blue	03020123
easyjet O-ring set	02900087
(1 small and 1 large O-ring for the spray head, 1 O-ring for the powder chamber cap)	



PERIO spray head



120° spray head



90° spray head



powder chamber cap

→ STARLIGHT UNO

Colorful and even easiest handling

LED curing light available in 6 different colours.
One touch handling, outstanding curing performance, exchangeable battery ... persuasive features, persuasive price!



→ STARLIGHT UNO

1 handpiece, 1 charging station

Color

- white
- grey
- turquoise
- yellow
- light blue
- violet

Ref. no.

- 05100092
- 05100093
- 05100094
- 05100095
- 05100096
- 05100097



→ TECHNICAL DATA

sizes (l x Ø in mm):

190 x 21

weight:

77 g

light source:

high power LED

wavelength:

dominant wavelength 440 – 465 nm

light intensity*:

> 1.500 mW/cm²

average life span LED:

1.800.000 cycles of 20 sec

battery:

rechargeable Lithium-ion battery
(no memory-effect)

capacity of the accumulator:

240 cycles of 10 sec

charging time:

4 hours (240 min)

warranty:

3 years, battery excluded

* measured with CURE RITE® radiometer by Dentsply

→ STARLIGHT PRO

Lightweight, powerful, versatile



→ STARLIGHT PRO

1 handpiece, 1 charging station

Color

white, shiny

yellow, shiny

silver

black, matt
carbon fiber look

Ref. no.

05100083

05100085

05100087

05100086

The lithium-ion battery of the cordless **starlight pro** provides 320 working cycles of 10 seconds – after only 120 minutes of charging. The power of more than 1.400 mW/cm² hardens a 2-mm thick layer in only 10 seconds. All this with only 105 g.



→ TECHNICAL DATA

sizes (l x Ø in mm):

190 x 23

weight:

105 g

light source:

high power LED

wavelength:

dominant wavelength 440 – 465 nm

optical fiber:

Ø 8 mm black, click-clack-connection, sterilisable at 135° C

light intensity*:

Ø 8 mm optical fiber, black > 1.400 mW/cm²

average life span LED:

1.800.000 cycles of 20 sec

accumulator:

lithium-ion-accumulator (no memory-effect)

capacity of the accumulator:

320 cycles of 10 sec

charging time:

about 120 min

warranty:

3 years (light guide and battery excluded)

* measured with CURE RITE® radiometer by Dentsply

→ TECHNICAL DATA

sizes (l x Ø in mm):

141 x 23

weight:

102 g

wavelength:

440 – 465 nm

optical fiber:

Ø 8 mm, sterilisable at 135° C

light intensity*:

Ø 8 mm optical fiber > 1.400 mW/cm²

average life span LED:

1.800.000 cycles of 10 sec

warranty:

3 years (light guide excluded)

→ STARLIGHT S⁺

Product

starlight s⁺

Ref. no.

05100098

* measured with CURE RITE® radiometer by Dentsply

→ STARLIGHT S⁺


flexibility integrated - 3 operating modes
easy handling

starlight s⁺ can be easily built into your dental unit and is simply placed in the handpiece holder – exactly where you need it.



STARLIGHT SX5

Built-in kit – the integrated solution

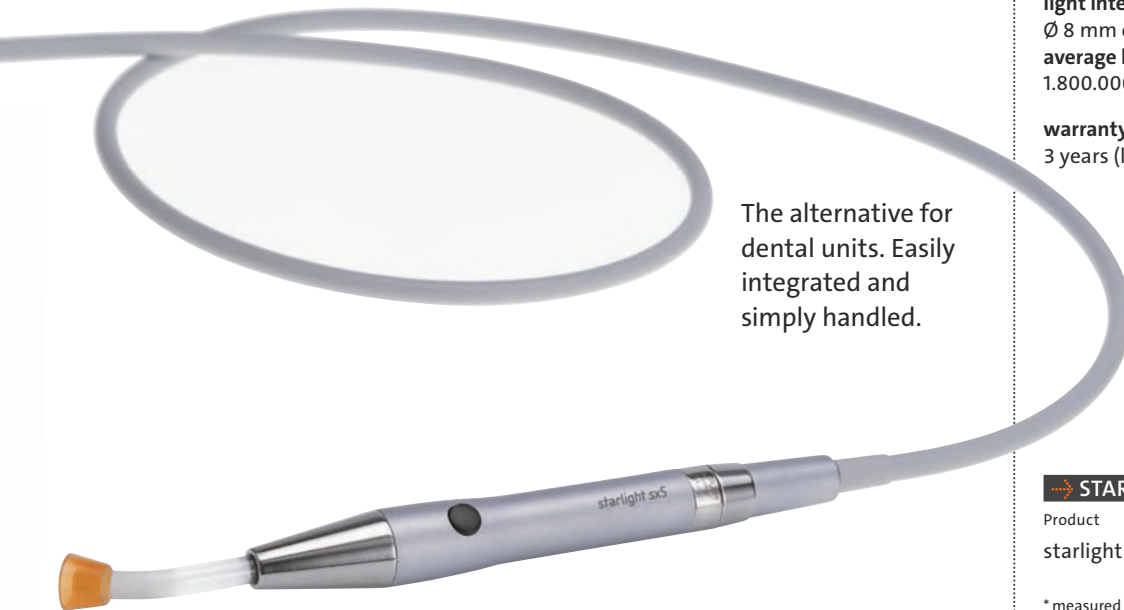

The alternative for dental units. Easily integrated and simply handled.

→ TECHNICAL DATA

sizes (l x Ø in mm):

148 x 22

weight

65 g

wavelength:

440 – 465 nm

optical fiber:

Ø 8 mm, click-clack-connection,
sterilizable at 135° C

light intensity*:

Ø 8 mm optical fiber > 1.400 mW/cm²

average life span LED:

1.800.000 cycles of 20 sec

warranty:

3 years (light guide excluded)

→ STARLIGHT SX5

Product

starlight sx5

Ref. no.

05100003

* measured with CURE RITE® radiometer by Dentsply



optical protection, 8mm



black optical fiber with metal click-clack connection Ø 8 mm



transparent optical fiber Ø 8 mm



black optical fiber Ø 8 mm



transparent optical fiber Ø 4,5 mm

SPARE PARTS

SPARE PARTS STARLIGHT UNO

Product

Ref. no.

battery starlight uno

03020222

protective sleeves (500 pieces)

03230015

optical protection (package with 2 pieces)

02900130

COMMON SPARE PARTS STARLIGHT PRO, S, P AND SLER AND BLUELIGHT

optical protection, sterilizable Ø 8 mm (package of 2 pcs.)

02900018

optical protection, sterilizable Ø 4.5 mm (package of 2 pcs.)

02900019

SPARE PARTS STARLIGHT PRO AND ORTHO with metal click-clack connection of the optical fiber

black optical fiber with metal click-clack connection Ø 8 mm

02900086

SPARE PARTS STARLIGHT PRO, S+ AND SX5

transparent optical fiber Ø 8 mm

02750009

black optical fiber Ø 8 mm

02750012

transparent optical fiber Ø 4.5 mm

02750013

starlight sx5 cord

02830029

SPARE PARTS STARLIGHT SLER®

transparent optical fiber with metal click-clack connection Ø 8 mm

02900085

black optical fiber with metal click-clack connection Ø 8mm

02900086



APPROPRIATE SETTINGS

INSERT	PIEZOSURGERY® <i>touch</i> PIEZOSURGERY® <i>white</i>	PIEZOSURGERY® 3	
INSERT	FUNCTION	MODE	POWER

IMPLANT PREP

IM1S	special*	BONE	SPECIAL
IM1 AL	implant	IMPL	IMPLANT
IM2A	implant	IMPL	IMPLANT
IM2A-15	implant	IMPL	IMPLANT
IM2P	implant	IMPL	IMPLANT
IM2P-15	implant	IMPL	IMPLANT
IM2.8A	implant	IMPL	IMPLANT
IM2.8P	implant	IMPL	IMPLANT
IM3A	implant	IMPL	IMPLANT
IM3A -15	implant	IMPL	IMPLANT
IM3P	implant	IMPL	IMPLANT
IM3P-15	implant	IMPL	IMPLANT
IM3.4A	implant	IMPL	IMPLANT
IM3.4P	implant	IMPL	IMPLANT
IM4A	implant	IMPL	IMPLANT
IM4P	implant	IMPL	IMPLANT
P2-3	implant	IMPL	IMPLANT
P3-4	implant	IMPL	IMPLANT

OSTEOTOMY

OT1	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT1A	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT2	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT4	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT5	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT5A	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT5B	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT6	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT7	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT7A	cortical/cancellous	BONE	CORTICAL/SPONGIOUS

INSERT	PIEZOSURGERY® <i>touch</i> PIEZOSURGERY® <i>white</i>	PIEZOSURGERY® 3	
INSERT	FUNCTION	MODE	POWER

OSTEOTOMY

OT7S-3	special*	BONE	SPECIAL*
OT7S-4	special*	BONE	SPECIAL*
OT7-20	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT8L	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT8R	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT9	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT11	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT12	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT12S	special*	BONE	SPECIAL*
OT13	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OT14	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
SLO-H	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
PL1	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
PL2	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
PL3	cortical/cancellous	BONE	CORTICAL/SPONGIOUS

OSTEOPLASTY

OP1	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP2	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP3	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP3A	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP4	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP5	perio	ROOT	PERIO
OP5A	perio	ROOT	PERIO
OP6	perio	ROOT	PERIO
OP7	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
OP8	special/perio	ROOT	PERIO
OP9	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
SLC	cortical/cancellous	BONE	CORTICAL/SPONGIOUS

FOR PIEZOSURGERY® INSERTS

INSERT	PIEZOSURGERY® <i>touch</i> PIEZOSURGERY® <i>white</i>	PIEZOSURGERY® 3	
INSERT	FUNCTION	MODE	POWER
SINUS MEMBRANE ELEVATION			
EL1	endo*	ROOT	ENDO*
EL2	endo*	ROOT	ENDO*
EL3	endo*	ROOT	ENDO*
SLS	special*	BONE	SPECIAL*
SLE1	special*	BONE	SPECIAL*
SLE2	special*	BONE	SPECIAL*
PERIOSTEUM PREPARATION			
PR1	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
PR2	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
ENDO RETRO			
EN1	endo*	ROOT	ENDO*
EN2	endo*	ROOT	ENDO*
EN3	endo*	ROOT	ENDO*
EN4	endo*	ROOT	ENDO*
EN5R	endo*	ROOT	ENDO*
EN5L	endo*	ROOT	ENDO*
EN6R	endo*	ROOT	ENDO*
EN6L	endo*	ROOT	ENDO*
EXTRACTION			
EX1	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
EX2	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
EX3	cortical/cancellous	BONE	CORTICAL/SPONGIOUS
EXL1	cortical	BONE	CORTICAL
EXL2	cortical	BONE	CORTICAL
EXL3	cortical	BONE	CORTICAL

INSERT	PIEZOSURGERY® <i>touch</i> PIEZOSURGERY® <i>white</i>	PIEZOSURGERY® 3	
INSERT	FUNCTION	MODE	POWER
EXPLANATION			
EXP3-R	special/cancellous/cortical	BONE	SPECIAL/SPONGIOUS/CORTICAL
EXP3-L	special/cancellous/cortical	BONE	SPECIAL/SPONGIOUS/CORTICAL
EXP4-R	special/cancellous/cortical	BONE	SPECIAL/SPONGIOUS/CORTICAL
EXP4-L	special/cancellous/cortical	BONE	SPECIAL/SPONGIOUS/CORTICAL
PERIO			
PP1	perio	ROOT	PERIO
PS1	perio	ROOT	PERIO
PS2	perio	ROOT	PERIO
PS6	perio	ROOT	PERIO
PP10	perio	ROOT	PERIO
PP11	perio	ROOT	PERIO
PP12	perio	ROOT	PERIO
IMPLANT CLEANING			
ICP + IC1	endo/perio	ROOT	ENDO/PERIO
RESTORATIVE			
DB2	special/perio	BONE ROOT	SPECIAL PERIO

* Maximum power permitted

→ APPROPRIATE SETTINGS



function	endo perio / scaler restorative	endo perio / scaler restorative	endo perio / scaler prosthetics	
----------	---------------------------------------	---------------------------------------	---------------------------------------	--

SCALING

S1-S	1-2**	1-2**	1-2**	0-30%
S1, S2, S3, S4, S5, S6	1-6**	1-5**	1-3**	0-100%

PERIO

P1, P2, P3, P4, P10, P11, P12, P13, P14	1-3**	1-3**	1-3**	0-50%
P15, P16R, P16L	1-3**	1-3**	1-3**	0-50%

IMPLANT CLEANING

ICS + IC1	1-2**	1-2**	1-2**	0-30%
-----------	-------	-------	-------	-------

ENDO

E1, E2	1-6	1-5	1-2**	0-10%
--------	-----	-----	-------	-------

FOR ULTRASOUND INSERTS



function	endo	perio/scaler	restorative	endo	perio/scaler	restorative	endo	perio/scaler	prosthetics
----------	------	--------------	-------------	------	--------------	-------------	------	--------------	-------------

ENDO REVISION

ER1		1-6		1-5			1-3		0-100%
ER2, ER3, ER4, ER5		1-3		1-3			1-3		0-50%

ENDO RETRO

RD3, RD4, RD5		1-3		1-3			1-3		0-50%
---------------	--	-----	--	-----	--	--	-----	--	-------

CROWN PREP

DB1			4-6		4-5			○**	50-100%
-----	--	--	-----	--	-----	--	--	-----	---------

EXTRACTION

ME1, ME2, ME3 ***			6		5				
-------------------	--	--	---	--	---	--	--	--	--

* These devices are not provided with a power level indicator. The percentage value indicated in the table refers to the adjustment percentage of the scale of the potentiometer supplied with the device.

** Can be used with or without "SOFT MODE"-function.

*** not available in the United States.


MECTRON ONLINE EDUCATION

We are happy to share with you our online education platform.

It will allow you to attend free of charge our live and on-demand educational webinars and view our clinical video library.

A fantastic possibility to take a deeper view on clinical opportunities helpful for your daily practice.

Visit our mectron education platform and register to receive notifications about upcoming webinars.

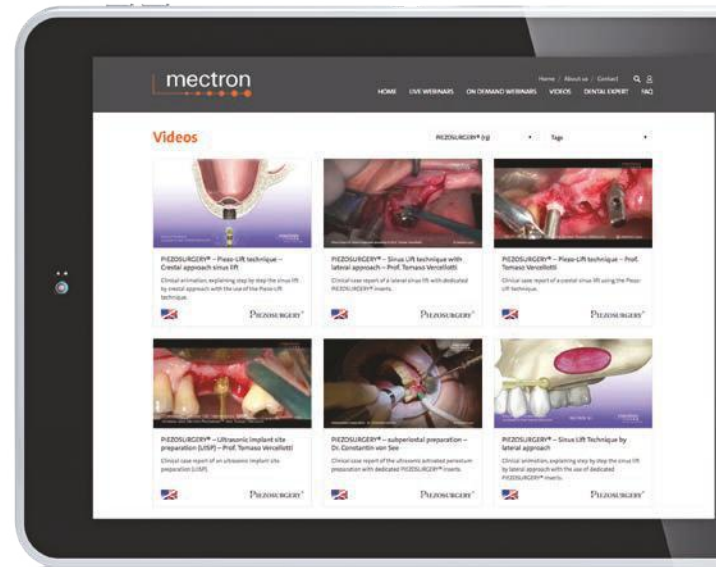



On Demand Webinars

- ON DEMAND WEBINAR** (UK)
Minimally invasive options for implant treatment of reduced bone ridges
Speaker: Dr. Claudio Stacci
REGISTER NOW
- ON DEMAND WEBINAR** (UK)
New implant shape to simplify narrow ridge treatment
Speaker: Prof. Soren Wenzel
REGISTER NOW
- ON DEMAND WEBINAR** (IT)
Evoluzione dei protocolli operativi della terapia parodontale non chirurgica: quali livelli di attenzione?
Speaker: Prof. Agg. Graziella Maria Nordi
REGISTER NOW WE PROPHYLAXIS
- ON DEMAND WEBINAR** (FR)
Gain d'attache parodontale, Gain en confiance des patients, sérénité retrouvée!
Speaker: Dr. Maria Eve Bezzone
REGISTER NOW WE PROPHYLAXIS
- ON DEMAND WEBINAR** (IT)
Opzioni mininvasive per il trattamento implantare delle creste ossee di spessore ridotto
Speaker: Dr. Claudio Stacci
REGISTER NOW
- ON DEMAND WEBINAR** (BR)
Aplicação clínica da piezocirurgia: benefícios e indicações
Speaker: Dr. Robert Carvalho da Silva
REGISTER NOW PILZBURGER

Live Webinars

- LIVE WEBINAR** (UK)
Oral skills in management: New technologies and methods
Speaker: Dr. Stefano Cecchi
REGISTER NOW WE PROPHYLAXIS
- LIVE WEBINAR** (DE)
Pubertäts Anwendung in der Prothetik – ist es Zeit?
Speaker: Sebastian Kapp
REGISTER NOW WE PROPHYLAXIS
- LIVE WEBINAR** (UK)
Contemporary aesthetic orthodontics
Speaker: Dr. Robert Carvalho da Silva
REGISTER NOW PILZBURGER



Free of charge registration
<https://education.mectron.com/>



DISTRIBUTORS

Please check also our permanently updated distributor list under www.mectron.com.

COUNTRY	DEALER	PHONE	E-MAIL	MECTRON PRODUCTS	PIEZOSURGERY® PRODUCTS
→ ALBANIA	Rodio Pharma	+355 68 2066757	info@rodiofarma.al	●	●
→ ALBANIA	Dental Leader Corsi	+355 69 2097616	mema_edi@yahoo.it	●	●
→ ALGERIA	Dentina	+213 02 3772515	dentinadental@hotmail.com	●	●
→ ARGENTINA	Biodental / Latinmarket SA	+54 11 43825012	jcdescals@latinmarket-sa.com.ar	●	●
→ AUSTRALIA	BioMeDent Pty Ltd	+61 1300 792624	info@biomedent.com.au	●	●
→ AUSTRIA	mectron Deutschland Vertriebs GmbH	+49 221 492015 0	info@mectron.de	●	●
→ AZERBAIJAN	ADENTAA Dental Supplier	+994 12 5251117	info@adentta.az	●	●
→ BELGIUM	Hofmeester Dental BV	+31 884 100318	info@robouw.nl	●	●
→ BOLIVIA	GENTE DENTAL	+32 67 844284	vegadent23@hotmail.com	●	●
→ BOSNIA-HERZEGOVINA	Mikodental BH d.o.o.	+381 15 345123	mikodental@gmail.com	●	●
→ BRAZIL	Implantec - Comércio e Assistência Técnica	+55 31 34983326	jarbas@implantecmg.com.br	●	●
→ BULGARIA	Medical Depot Rik Ood	+359 2 8052500	office@mdrk.bg	●	●
→ CANADA	Piezosurgery® Inc. / Mectron US	+1 614 459 4922	info@piezosurgery.us	●	●
→ CHILE	Berghem	+56 2 2144643	info@berghem.cl	●	●
→ CHINA	Shanghai Tian Ying Medical Instrument Co., Ltd.	+86 21 51060620	sales@tydental.com	●	●
→ COLOMBIA	REPRESENTACIONES A&D LTDA.	+57 9156636	info@ayd.com.co	●	●
→ COSTA RICA	ODONTOWEB PROFESIONAL SA	+506 22 34 31 33	info@odontoweb.com	●	●
→ CROATIA	Dental Grupa D.O.O.	+385 52 752106	info@dentalgrupa.hr	●	●
→ CZECH REPUBLIC	mectron.cz s.r.o.	+420 777 761107	info@mectron.cz	●	●
→ DENMARK	Eltident APS	+45 70271747	tina@eltident.com	●	●
→ ECUADOR	PRODENTEC CIA. LTDA.	+593 72831232	ventas@prodentec.com.ec	●	●
→ EGYPT	Zaident for import dental supply	+20 111 222 5076	zaident@link.net	●	●

COUNTRY	DEALER	PHONE	E-MAIL	MECTRON PRODUCTS	PIEZOSURGERY® PRODUCTS
→ ESTONIA	Prodent OÜ	+372 656 5982	info@prodent.ee	●	●
→ ESTONIA	Retent	+372 7 301 400	retent@retent.ee	●	●
→ FINLAND	Plandent Oy	+358 20 77952210	markku.nappari@plandent.com	●	●
→ FRANCE	mectron France	+39 0185 1759069	mectronfrance@mectron.fr	●	●
→ GEORGIA	UNIMEDI LTD	+995 595900344	info@unimedi.ge	●	●
→ GERMANY	mectron Deutschland Vertriebs GmbH	+49 221 492015 0	info@mectron.de	●	●
→ GREECE	ARISTON DENTAL	+30 694 8622608	info@aristondental.com	●	●
→ GREECE	Green Dental lke	+30 210 3803055	siakavaras.g@gmail.com	●	●
→ GUATEMALA	IMFOHSA	+502 2317 2100	info@imfohsa.com	●	●
→ HONG KONG	Denco (HK) Medical Supplies Ltd.	+852 2493 7112	hugo@denco.com.hk	●	●
→ HUNGARY	Dental Unio Kft	+36 96 518047	info@dentalunio.hu	●	●
→ HUNGARY	Medi-cont Kft.	+36 96 516360	info@medi-cont.hu	●	●
→ INDONESIA	PT. Andini Sarana	+62 21 4610515	hadiwira@andini.co.id	●	●
→ IRAN	Nik Nam Teb Ayandeh	+98 21 66467007	info@niknamteb.com	●	●
→ IRELAND	General Medical	+44 1380 734990	info@generalmedical.co.uk	●	●
→ ISRAEL	Bicon Isreal	+972 86450479	biconisrael@gmail.com	●	●
→ ISRAEL	E. Dental Ltd.	+972 36915507	elih@e-dental.co.il	●	●
→ ITALY	mectron s.p.a.	+39 0185 35361	mectron@mectron.com	●	●
→ JAPAN	Implatex Co., Ltd	+81 3 58508555	itx@itx.co.jp	●	●
→ JAPAN	Tokyo Dental Industrial Co. Ltd.	+81 3 38310176	orii@tokyodental.co.jp	●	●
→ JORDAN	Ibn Rushd Dental Supplies (IRDS)	+962 655 26162	tareq@irds.jo	●	●
→ KUWAIT	Al-Sayafe Medical & Pharmacy Supplies	+965 2245 4924	alsayafe@alsayafemedical.com	●	●



DISTRIBUTORS

COUNTRY	DEALER	PHONE	E-MAIL	MECTRON PRODUCTS	PIEZOSURGERY® PRODUCTS
→ LATVIA	AS "Recipe Plus"	+371 678 1584242	info@recipe.lv	●	
→ LEBANON	Interbio SAL	+961 1 491891	chadi.maroun@interbio.me	●	●
→ LITHUANIA	Jsc Med Grupe	+370 5 2761377	giedre@medgrupe.lt	●	●
→ LITHUANIA	Vilimekso Servisas JSC	+370 5 2796709	info@vilimeksas.lt	●	
→ LUXEMBOURG	T.D.L.s.a.	+352 446881	info@tdlsa.net	●	●
→ MALAYSIA	Galla Dental Cooperation Sdn Bhd	+603 6259 3800	info@galladental.com		●
→ MEXICO	DENTADEC S.A. DE C.V.	+55 55454520	daniel@dentadec.com		●
→ MOLDOVA	VIVODENT SRL CENTRUL STOMATOLOGIC	+373 22 226522	office@vivodent.md	●	●
→ MONTENEGRO	MIKODENTAL d.o.o.	+381 15 345123	mikodental@gmail.com	●	●
→ MOROCCO	Dental Project	+212 522 47 37 13	Sales@dentalproject.ma	●	●
→ NETHERLANDS	Hofmeester Dental Bv	+31 884 100318	verkoop@hofmeester.nl	●	●
→ NICARAGUA	INDENT S.A.	+505 2223 3517	info@indentsa.com	●	●
→ NORWAY	Dentalvarehuset Spesialistsenter AS	+47 22 333900	thomas@dentalvarehuset.no	●	●
→ NORWAY	Technomedics Norge AS	+47 6988 7920	ek@technomedics.no	●	●
→ PAKISTAN	Glow Pak International	+92 42 58460902	info@glowpak.com	●	●
→ POLAND	FM Produkty Dla Stomatologii	+48 12 632 19 27	biuro@fmdental.pl	●	●
→ PORTUGAL	Incotrading S.A.	+34 91 3807490	incotrading@incotrading.net	●	●
→ ROMANIA	S.I.T.E.A. Romania Srl	+40 215 280320	fabrizio@sitea.ro	●	●
→ RUSSIA	Anas Ltd	+7 495 626 5455	clients@anas.ru	●	●
→ SAUDI ARABIA	Cigalah healthcare company	+966 12 231 6194	ifarra@cigalah.com.sa	●	●
→ SERBIA	MIKODENTAL Depot	+38115 345123	mikodental@gmail.com	●	●
→ SINGAPORE	SA Dental Supply PTE LTD	+65 6295 2568	info@sadental.com.sg	●	●

COUNTRY	DEALER	PHONE	E-MAIL	MECTRON PRODUCTS	PIEZOSURGERY® PRODUCTS
→ SLOVAKIA	mectron.cz s.r.o.	+420 777 761107	info@mectron.cz	●	●
→ SLOVENIA	Prodent Int. D.O.O.	+386 1 5204800	info@prodent.si	●	●
→ SOUTH KOREA	MD Safe Co. Ltd.	+82 2 334 2815	jclosa@naver.com	●	●
→ SPAIN	Incotrading S.A.	+34 91 3807490	incotrading@incotrading.net	●	●
→ SWEDEN	Dentalmind	+46 346 48800	marcus@dentalmind.com	●	●
→ SWITZERLAND	Bien-Air Dental SA	+41 032 344 64 64	sales-ch@bienair.com	●	●
→ SYRIA	EBDA Engineering & Development	+963 1157 34444	management@ebda-eng.net	●	●
→ TAIWAN	Sudent Corporation	+886 2822 80 861	sudent.co@msa.hinet.net	●	●
→ THAILAND	Udom Medical Equipment Co. Ltd.	+662 322 7124	wisit@udom.co.th	●	●
→ TUNISIA	New Med Plus	+216 31 561 563	newmedplus@gmail.com	●	●
→ TURKEY	Tekmil Tibbi Arac	+90 21 2337837	info@teknil.com.tr	●	●
→ UKRAINE	Galit	+380 352 433807	office@galit.te.ua	●	●
→ UNITED ARAB EMIRATES	New Al Farwaniya Surgicals & Medical Equipment L.L.C.	+971 2 6775447	info@alfarwaniya.com	●	●
→ UNITED KINGDOM	General Medical	+44 1380 734990	info@generalmedical.co.uk	●	●
→ USA	Piezosurgery® Inc./Mectron US	+1 614 459 4922	info@piezosurgery.us	●	●
→ VIETNAM	Pacific Dental Supply	+84 8 391121270	customer-care-hcm@dentalpacific.com.vn	●	●

mectron s.p.a.,
via Loreto 15/A, 16042 Carasco (Ge), Italia,
tel +39 0185 35361, fax +39 0185 351374

→ www.mectron.com or mectron@mectron.com



© Copyright mectron S.p.A., Carasco, Italy
All rights reserved. Texts, pictures and graphics of mectron brochures are protected by copyright and other protection laws. Without written approval of mectron S.p.A. the contents may not be copied, distributed, changed or made available to third parties for commercial purposes.

Please consult product labels and manual for indications, contraindications, hazards, warnings, cautions.
Please contact your local sales representative for product availability.

The following products are registered under MDR 2017/745. For graphic reasons the styling of the product names might be different. The registered product names correspond as following:

- | | | | |
|-----------------------------|----------------------------|--------------------------|------------------------|
| • PIEZOSURGERY <i>touch</i> | PIEZOSURGERY® <i>touch</i> | • COMPACT PIEZO LED | COMPACT PIEZO LED |
| • PIEZOSURGERY <i>touch</i> | PIEZOSURGERY® <i>touch</i> | • COMPACT PIEZO P2K | COMPACT PIEZO P2K |
| • combi <i>touch</i> | COMBI <i>touch</i> | • starjet | STARJET |
| • multipiezo | MULTIPIEZO | • SUBGINGIVAL PERIO TIPS | Subgingival Perio Tips |
| • MICROPIEZO | MICROPIEZO | | |